


DuraSeal® Dural Sealant System

During neurosurgical procedures, it is essential to prevent the leakage of cerebrospinal fluid (CSF) in order to reduce potential adverse events and ensure positive patient outcomes. DuraSeal® dural sealant system, hydrogel made with a specific polymer technology, provides a quick and effective watertight seal intraoperatively and through the critical healing period

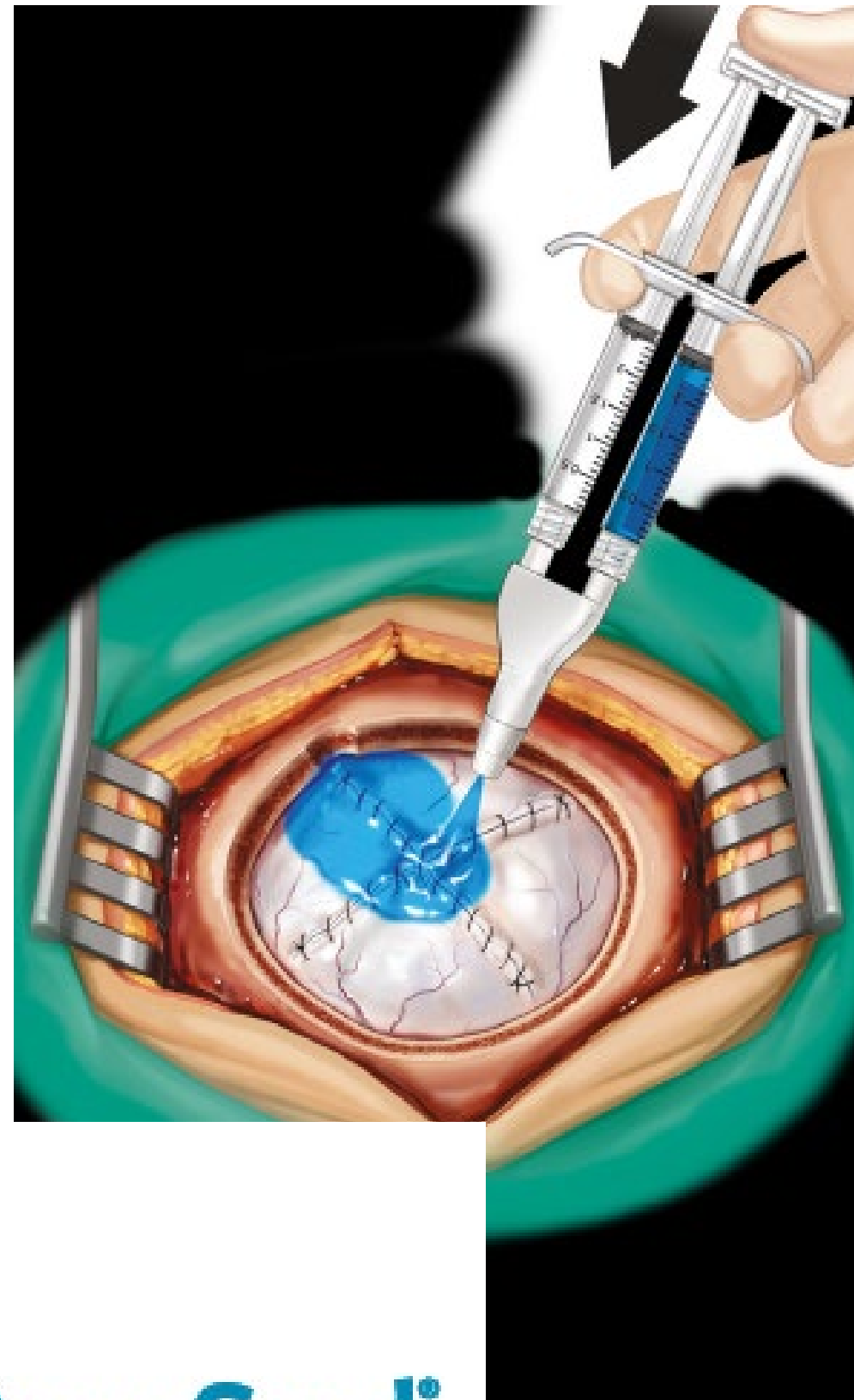

STRENGTH

Tissue adherence and cohesive strength to withstand critical pressures.

BIOCOMPATIBILITY

Biocompatible PEG (Polyethylene glycol) hydrogel.

VISIBILITY


Distinctive blue colorant provides visualization to assess sealant coverage and thickness.

SPEED

Prepared in less than two minutes. When applied, forms a watertight seal in seconds.

CONVENIENCE


Single kit, stored at room temperature.

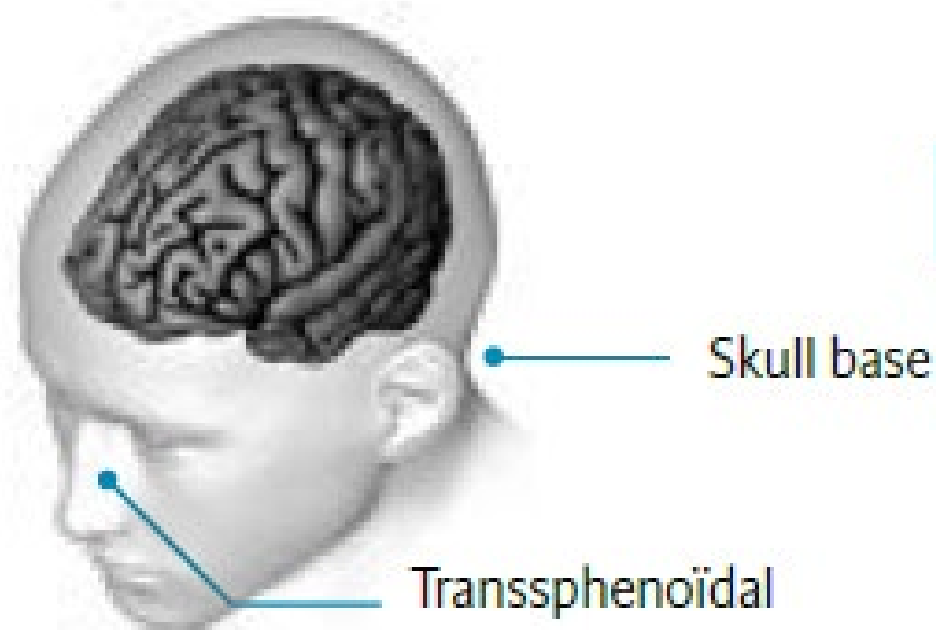

DuraSeal®
Cranial Sealant System

DuraSeal® Xact spinal sealant system to ensure a watertight dural closure

During spinal surgery procedures, it is essential to prevent the leakage of cerebrospinal fluid (CSF) in order to reduce potential adverse events and ensure positive patient outcome. DuraSeal® Xact sealant system provides surgeons with a polymer technology that provides a watertight seal quickly and effectively, intraoperatively and through the critical healing period.

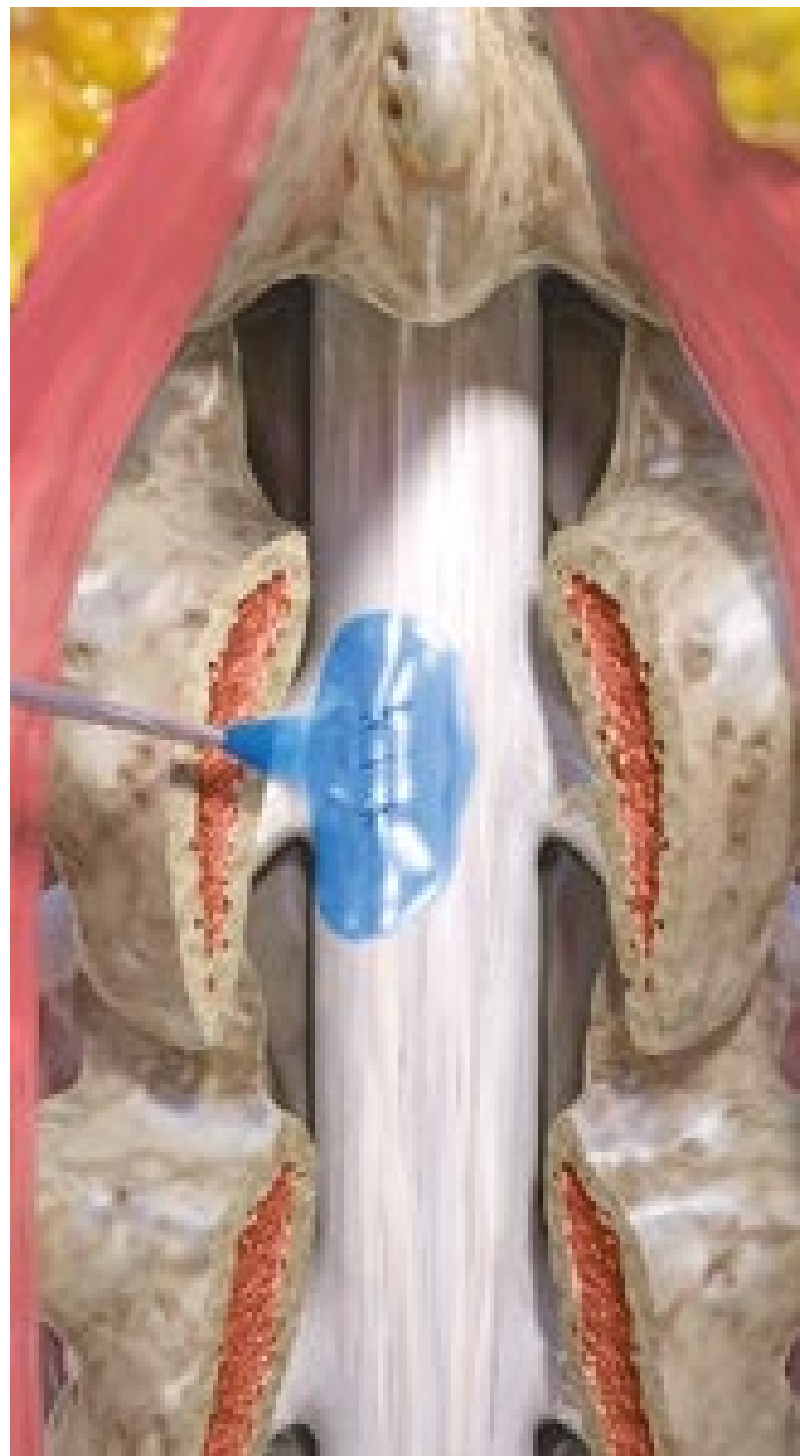

STRENGTH

Tissue adherence and cohesive strength to withstand critical pressures.

BIOCOMPATIBILITY

Biocompatible PEG (Polyethylene glycol) hydrogel.

VISIBILITY

Distinctive blue colorant provides visualization to assess sealant coverage and thickness.

SPEED

Prepared in less than two minutes. When applied, forms a watertight seal in seconds.

FEATURE

Limited swelling formula mitigates concern of expansion in confined spaces.

Fully synthetic and absorbable in 9-12 weeks for normal healing.

CONVENIENCE

Single kit, stored at room temperature.

Extended Tip Applicator

The Extended Tip Applicators give DuraSeal® Xact sealant system users the versatility of a malleable manual applicator with extended reach and added visibility. Available in 8 cm and 15 cm length.

60° bendable shaft shown with 8 cm Extended Tip Applicator, 205108



- Enhance your approach in hard-to-reach areas of the spine and brain
- Offer flexibility that extends access, through your surgical approach, with a 60° bendable shaft
- Pause and resume application of DuraSeal® and DuraSeal® Xact cranial and spinal sealant systems intra-operatively

MicroMyst® Applicator

The flexible air-assisted MicroMyst® applicator delivers precise application through a fine mist spray (14 cm length). Used with the Flow Regulator for the controlled application of two liquids.

